

3 1761 11556784 4




CAI HW190

58 H55



HOME NURSING TRAINING KIT



Digitized by the Internet Archive
in 2022 with funding from
University of Toronto

<https://archive.org/details/31761115567844>

CAI HW 190
58H55

HOME NURSING TRAINING KIT

**Produced for Civil Defence
by Information Services Division**

**Department of National Health
and Welfare**

**by authority of the Minister
Hon. J. Waldo Monteith**

1958

HOME NURSING TRAINING KIT

BOXES 1 AND 2


This illustrated folder lists the contents of the two boxes which make up a Home Nursing Training Kit, and gives explicit instructions regarding repacking and shipping. These instructions were written to aid the lecturer and students as to the upkeep of equipment and supplies and the replacement of the items in the correct apertures, thus ensuring the maximum life to and the proper care of the equipment during use and shipping. Numbers on the illustrations indicate the correct apertures for all equipment.



The contents of this Kit should be kept clean at all times. It would be appreciated if the instructors of Home Nursing classes make arrangements for the laundering of all used linen prior to the shipment of the Kit to a forwarding address.

SHIPPING INSTRUCTIONS

1. The locks are universal and may be opened with any one of the keys shipped with the Home Nursing Training Kits. The keys are to be shipped under separate cover to the specific address of the place of instruction as supplied by the teaching organization.
2. All shipping instructions and the addresses are to be placed on a shipping tag attached to the hasp on each box. The outside of the boxes are not to be defaced with shipping instructions or by any other information.
3. This Kit is the property of the Civil Defence Health Services and is under the control of the Provincial Civil Defence Co-ordinator who will be responsible for the payment of all shipping charges. The Kits are to be forwarded prepaid by the Provincial Civil Defence Co-ordinator to the addresses at which the instruction will take place, which will be supplied by the teaching organization. They are to be returned by the teaching organization, collect, to the address indicated by the Provincial Civil Defence Co-ordinator.




<i>Aperture No.</i>	<i>Description</i>	<i>Quantity</i>	<i>Instructions</i>
1	Bed, single, roll away type	1	The mattress is to be covered at all times by the mattress cover provided. The bed and mattress are to be assembled as illustrated before returning them to Box No. 1.
	Mattress	1	
	Mattress Cover	1	

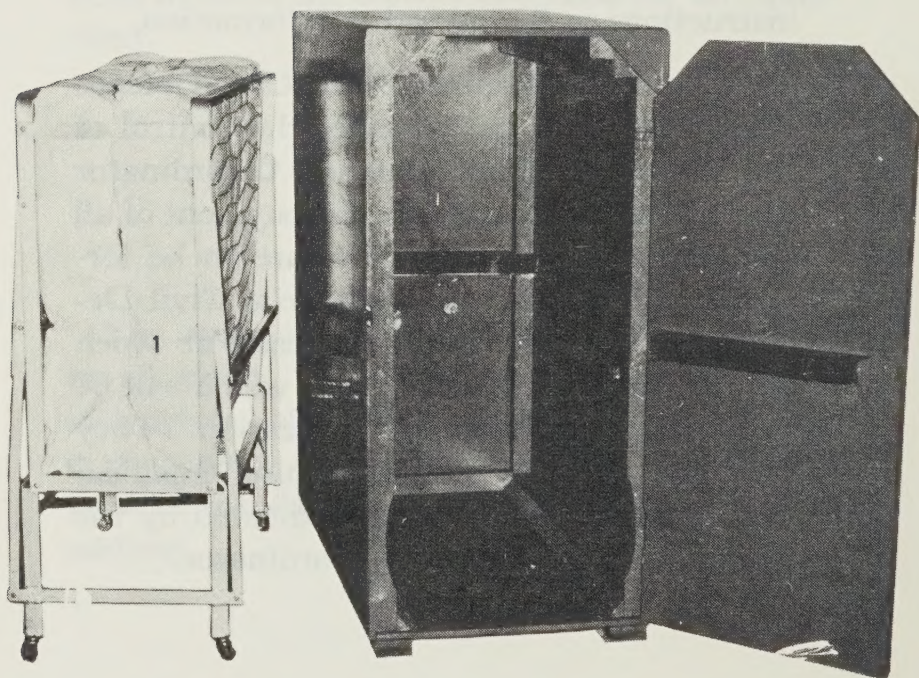

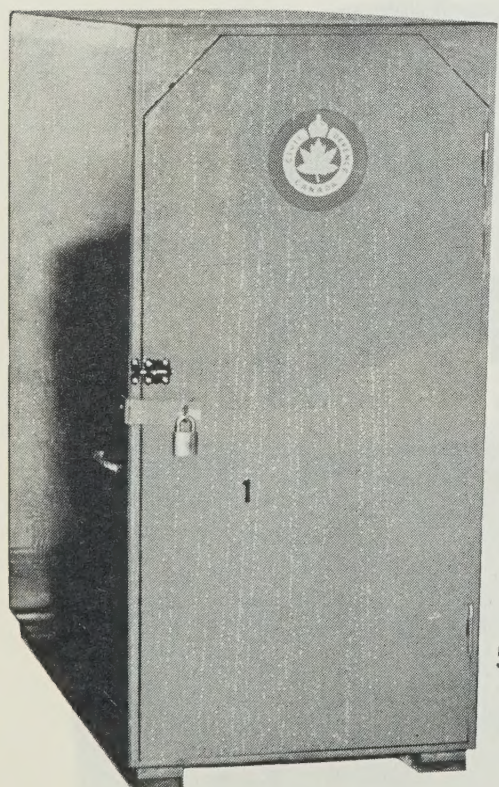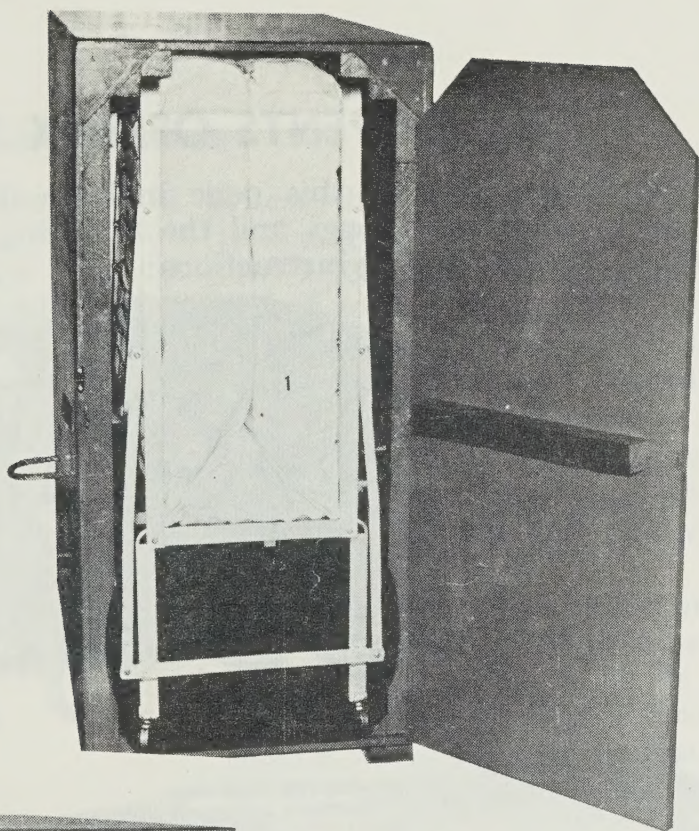




CONTENTS OF BOX 2

The illustrations on this page indicate the apertures to which each item belongs and the following information is more specific packing instructions:

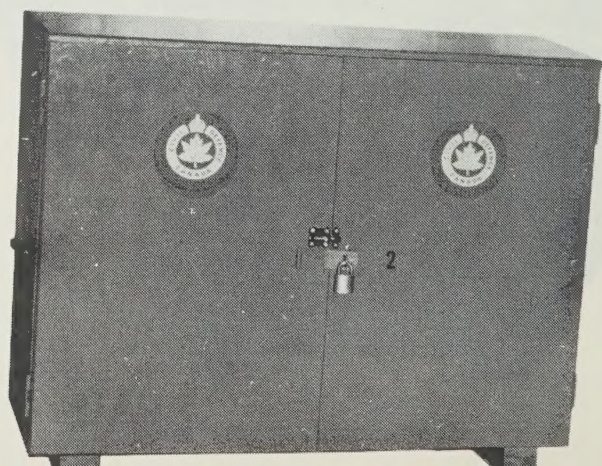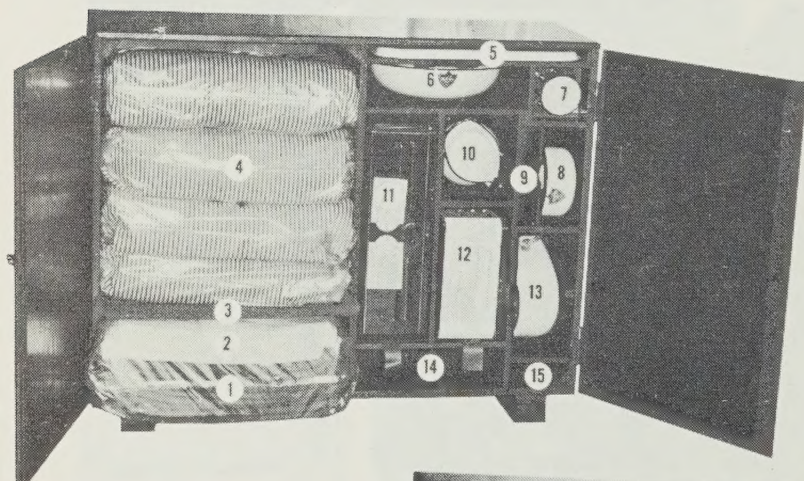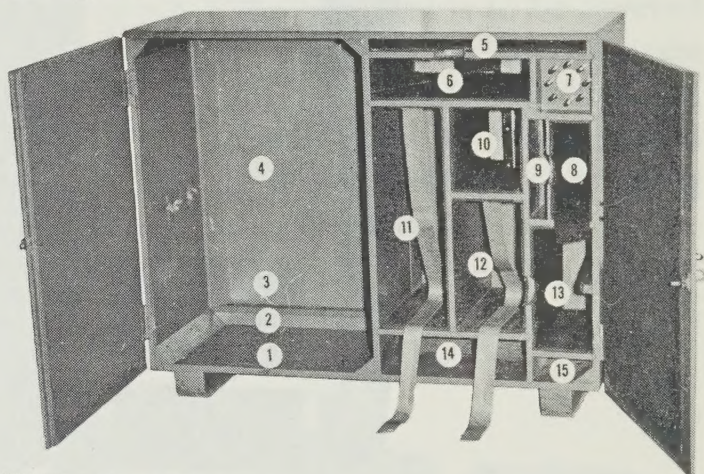

<i>Aperture No.</i>	<i>Description</i>	<i>Quantity</i>	<i>Instructions</i>
1	Blankets, Grey	2	{ These three items are to be placed in a polyethylene bag (36" x 22"), the blankets put in first so that they would be on the bottom. The contents of this bag are then placed in Box No. 2 as illustrated.
2	Sheets, white	7	
	Pillow Cases, White	4	
3	Rubber Sheeting	1	This item is unwrapped and is to be placed in Box No. 2 on top of the sheets and pillow cases as illustrated.
4	Pillows	4	The pillows are to be placed in the individual polyethylene bag and the air is to be deflated from the bag before closing the end of it. The pillows are then to be placed in Box No.2 as illustrated.
5	Tray, Serving, W.E.	1	This item must be inserted upside down so that it may be easily withdrawn by means of the cutout portion of the shelf.
6	Basin, Round Wash, W.E.	1	This item must be inserted right side up so that advantage may be taken of the cutout portion of the shelf.
7	Thermometers, oral, in bakelite case.	6	The thermometers are to be carefully returned to the bakelite cases. These cases are to be fitted into the holes arranged in a circular manner.
	Thermometers, rec- tal, in bakelite case.		
	Funnel, W.E.,	1	
	Catheter tip 2 oz.		The tip of the funnel is to be fitted into the centre hole. (See illustrations).

CONTENTS OF BOX 2 (Continued)

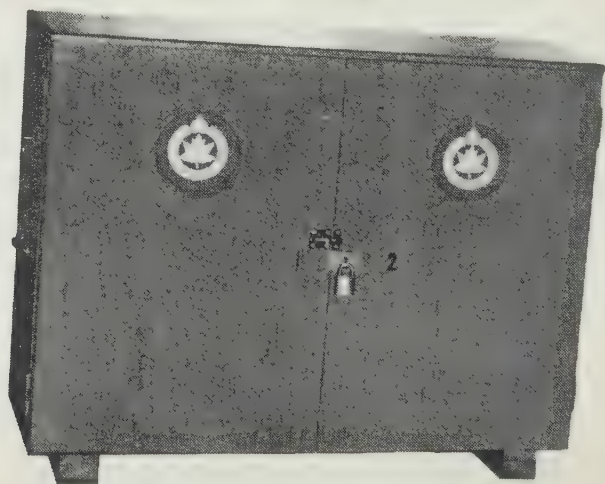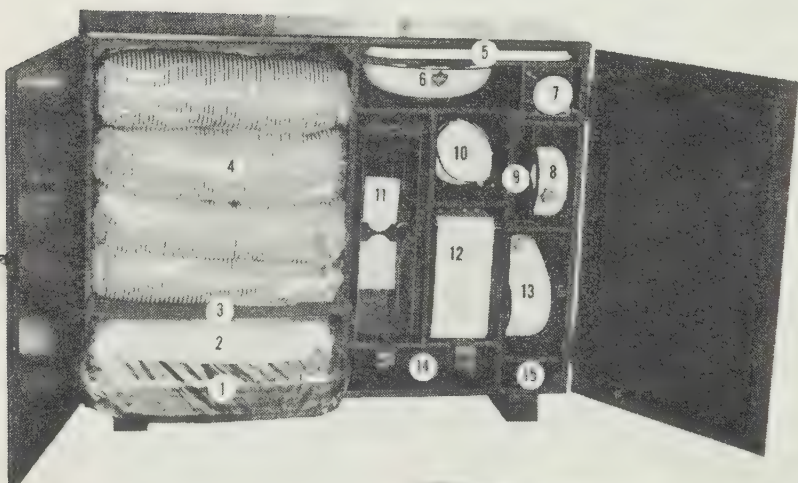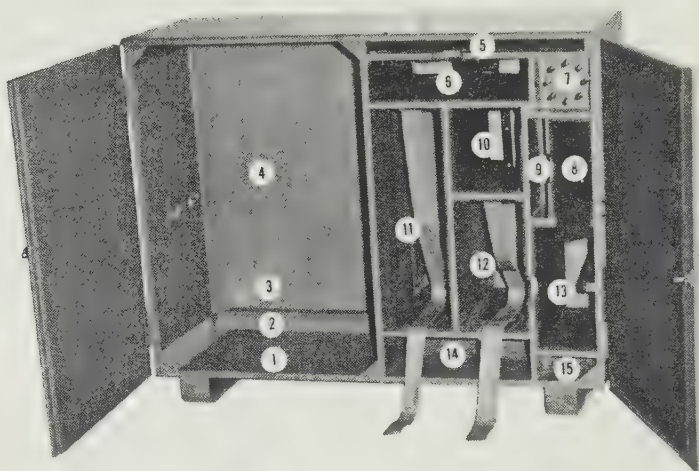

<i>Aperture No.</i>	<i>Description</i>	<i>Quantity</i>	<i>Instructions</i>
8	Basins, round, W.E.	2	The first basin is to be inserted, then the felt tongue is put into place. When the second basin is placed into the section, this leaves the felt tongue between the two basins. The two basins should be inserted so that the tops of the basins are next to the cut-out portion of the divider, as illustrated.
9	Instruction Booklet	1	This section has been reserved for the instruction booklet which is to be returned with the unit, as illustrated.
10	Pitcher, water, with handle, W.E. size 3 qts.	1	This item must be inserted with the open end out so that the handle will be readily available to assist with the removing of it, as illustrated.
11	Hot plate, two burner	1	The removing straps for these two items are clearly illustrated. These straps must be
12	Hot plate, one burner	1	behind the hot plates when they are returned to the apertures indicated for them. A small portion of each strap must remain at the bottom so that it may be used to remove the hot plates from the box.

CONTENTS OF BOX 2 (Concluded)

<i>Aperture No.</i>	<i>Description</i>	<i>Quantity</i>	<i>Instructions</i>
13	Bed Pan, W.E.	1	The hot water bottle; kidney basin; attachment set and rectal tube are to be wrapped and placed inside of the bed pan. Then the bed pan is to be inserted into the section with the top facing the cut-out portion of the divider, as illustrated.
	Bottle, hot water size 3 qt.	1	
	Set, attachment complete	1	
	Tube, rectal, red rubber, size 28, Fr.	1	
	Basin, Dressing, W.E.K.—shaped 8"	1	
13	Clamp, enemata	1	This item should be wrapped in tissue paper and stored with the attachment set in the bed pan.
14	Blank Section		
15	Blank Section		
4a	Gown white	1	The five items are to be placed in a polythene bag, the bath blanket first so that it will be on the bottom. The contents of this bag are then placed above aperture No. 4. (see illustrations)
	Washcloths	2	
	Towel, bath	1	
	Towel, face	1	
	Blanket, bath (flannelette sheet)	1	

ADDITIONAL EQUIPMENT

BOX No. 1

<i>Aperture No.</i>	<i>Description</i>	<i>Quantity</i>	<i>Instructions</i>
1	Bed Blocks (7"x8"x6")	4	These blocks are to be used to raise the bed to a more workable height. The smaller ones should be placed under the castors to raise the centre of the bed, and the larger ones are to be used to raise the ends of the bed. Store in space under bed in Aperture No. 1.
1	Bed Blocks (7"x5"x6")	4	

NOTE:

The supplies contained in these two boxes are the more expensive basic items and those most difficult to assemble for home nursing teaching. On the following page we have listed, for the information of instructors, additional articles which would have to be procured locally.

ARTICLES TO BE OBTAINED LOCALLY

Soap	Infant feeding bottle and nipple
Paper towels	Bottle for specimen or sputum
Wax paper	Absorbent
Newspapers	Adhesive
Toilet paper	Gauze and old cotton (for mustard paste)
Paper tissues	Dishes for invalid tray
Powder	Knife
Rubbing alcohol	Large spoon
Mustard	Large baby doll for infant demonstrations, nightgown
Flour	Diapers
Salt	Safety pins
Soda bi-carbonate	Hair comb, fine comb
Medicine bottle (filled with coloured water)	Tooth brush and tumbler
Medicine glass	Minute hand watch
Medicine dropper	Drinking tube or straws
Wooden applicator	Water pail and slop pail (if no running water is available)
Spring clothes-pins or forceps	Blackboard
Box of pills (soda bi-carbonate)	

NOTE:

In the most recent issue of home nursing kits the roll-away bed has been adjusted to a workable height, therefore, no bed blocks are included in this new issue. Otherwise the packaging of this new issue is exactly the same as previously described in this brochure.

THE QUEEN'S PRINTER AND CONTROLLER OF STATISTICS
OTTAWA, 1958

Cat. No. H84-1758

Ek 26.4.61

**University of Toronto
Library**

**DO NOT
REMOVE
THE
CARD
FROM
THIS
POCKET**

**Acme Library Card Pocket
LOWE-MARTIN CO. LIMITED**

